
### Initial Interrogation: Quick Look II

Device: Advisa DR MRI A3DR01 Serial Number: PZK875895S Date of Visit: 07-Mar-2018 13:02:36 Patient: Tsou, Chao Tsung/HA59891 ID: J120445951 Physician: Liao Min Tsun 0939429152 History: Sinus Node Dysfunction, Occasional 2° or 3° AV Block Device Status (Implanted: 08-Dec-2017) Battery Voltage (RRT=2.83V) (07-Mar-2018) 3.03 V Remaining Longevity 9 years (7.5 - 10 years) (based on initial interrogation) Atrial(5076) RV(5076) Lead Impedance 570 ohms 570 ohms Capture Threshold 0.625 V @ 0.40 ms 0.750 V @ 0.40 ms Measured On 07-Mar-2018 07-Mar-2018 2.00 V / 0.40 ms Programmed Amplitude/Pulse Width 1.50 V / 0.40 ms Measured P/R Wave 1.5 mV 12.3 mV Programmed Sensitivity 0.30 mV 0.90 mV **Parameter Summary** Mode DDD Lower Rate 60 bpm Paced AV 180 ms Mode Switch 171 bpm Upper Track 130 bpm Sensed AV 150 ms Upper Sensor 130 bpm Detection Rates **Therapies** AT/AF Monitor >171 bpm All Rx Off VT Monitor >150 bpm **Clinical Status** Since 08-Dec-2017 Longest Cardiac Compass Trends (Dec-2017 to Mar-2018) Treated AT/AF(Monitor) AT/AF 60 Monitored VT (>4 beats) 0 Fast A&V 0 20-AT/AF 2 61 sec 0 Time in AT/AF <0.1 hr/day (<0.1%) >8 Patient Activity **Functional Last Week** 6-(hr/day) Patient Activity 4.1 hr/day 4-2 0. Jan-18 Mar-18 May-18 Jul-18 Sep-18 Nov-18 Jan-19 Pacing (% of Time Since 08-Dec-2017) AT/AF Therapy Summary Pace-Terminated Episodes 0 Total VP 99.5 % (MVP Off) AS-VS 0.5 % AS-VP 89.8 % AP-VS <0.1 % AP-VP 9.7 %

### OBSERVATIONS (0)


<sup>·</sup> No observations based on current interrogation.

Initial Interrogation: Cardiac Compass Trends
Serial Number: PZK875895S Device: Advisa DR MRI A3DR01 Date of Visit: 07-Mar-2018 13:02:36 ID: **J120445951** Patient: Tsou, Chao Tsung/HA59891 Physician: Liao Min Tsun 0939429152

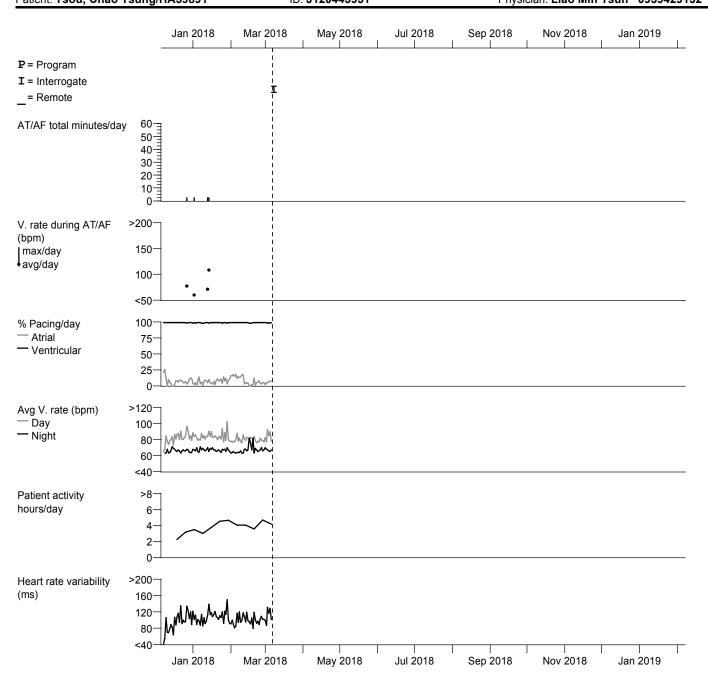

Initial Interrogation: Cardiac Compass Trends
Serial Number: PZK875895S Device: Advisa DR MRI A3DR01 Date of Visit: 07-Mar-2018 13:02:36 ID: **J120445951** Patient: Tsou, Chao Tsung/HA59891 Physician: Liao Min Tsun 0939429152

OptiVol fluid index is an accumulation of the difference between the daily and reference impedance.

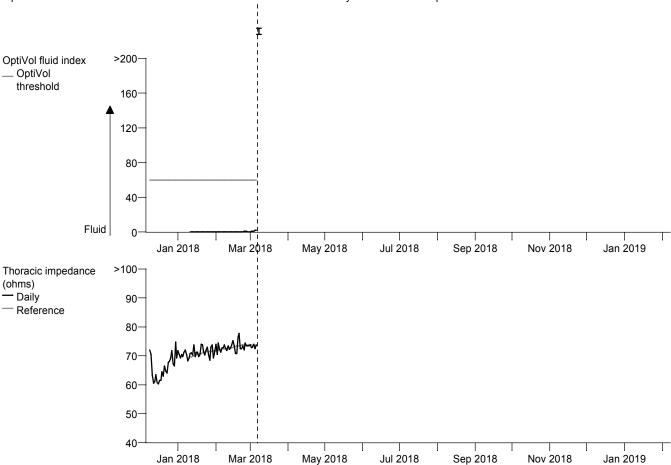

Serial Number: PZK875895S Device: Advisa DR MRI A3DR01 Date of Visit: 07-Mar-2018 13:02:36 ID: **J120445951** Patient: Tsou, Chao Tsung/HA59891 Physician: Liao Min Tsun 0939429152

#### **Prior to Last Session Since Last Session** No Prior Session 08-Dec-2017 to 07-Mar-2018 89 days % of Time **Total VP** 99.1 % AS-VS 0.5 % AS-VP 89.8 % AP-VS < 0.1% AP-VP 9.7 % Atrial % of Time 100-100-80 80 $\square$ AS 60-60-No Data Available ■AP 40 40 20 20 0 60 80 100 120 140 160 180 200 220> 60 80 100 120 140 160 180 200 220> < 40 Atrial Rate (bpm) Atrial Rate (bpm) Ventricular 100-100-% of Time 80-80-□VS 60-60-No Data Available ■VP 40-40-20 20 0-0. 60 80 100 120 140 160 180 200 220> < 40 60 80 100 120 140 160 180 200 220> Ventricular Rate (bpm) Ventricular Rate (bpm) Ventricular Time in AT/AF = 71 seconds 100-100-Rate **During** 80-80 AT/AF 60-60-No Data Available % of AT/AF 40-40 □VS 20 20 ■VP 80 100 120 140 160 180 200 220> 80 100 120 140 160 180 200 220> 60 Ventricular Rate (bpm) Ventricular Rate (bpm)

## Initial Interrogation: VT/VF Counters Serial Number: PZK875895S

Device: Advisa DR MRI A3DR01 Serial Number: PZK875895S Date of Visit: 07-Mar-2018 13:02:36
Patient: Tsou, Chao Tsung/HA59891 ID: J120445951 Physician: Liao Min Tsun 0939429152

VT/VF Counters	<b>Prior Session</b> No Prior Session	08-De	Session ec-2017 to ar-2018 lys
VT (>150 bpm)		0	
VT-NS (>4 beats)		0	
Fast A&V		0	
PVC Runs (2-4 beats)		<0.1	per hour
PVC Singles		16.2	per hour
Runs of VRS Paces		0.0	per hour
Single VRS Paces		0.0	per hour

Initial Interrogation: VT/VF Counters
Printed: 07-Mar-2018 13:23:22

Device: Advisa DR MRI A3DR01 Serial Number: PZK875895S Date of Visit: 07-Mar-2018 13:02:36
Patient: Tsou, Chao Tsung/HA59891 ID: J120445951 Physician: Liao Min Tsun 0939429152

AT/AF Summary	<b>Prior Session</b> No Prior Session	08-0	t Session Dec-2017 to Mar-2018 ays
% of Time AT/AF		<0.1	%
Average AT/AF time/day		<0.1	hours/day
Monitored AT/AF Episodes		<0.1	per day
Treated AT/AF Episodes		0.0	per day
Pace-Terminated Episodes		0.0	%
% of Time Atrial Pacing		9.7	%
% of Time Atrial Intervention		0.0	%
AT-NS (>6 beats)		<0.1	per day

### Since Last Session 08-Dec-2017 to 07-Mar-2018

AT/AF Du	rations				AT/	AF Episo	de Start 1	Γime		
Duration	Episodes	% Epis	odes							
>72 hr	0	100								
48 hr to 72 hr	0									
24 hr to 48 hr	0	80-								
12 hr to 24 hr	0									
4 hr to 12 hr	0	60-								
1 hr to 4 hr	0									
10 min to 1 hr	0	40-								
1 min to 10 min	1									
<1 min	4	20-								
		0+	09:00	12:00	15:00	18:00	21:00	00:00	03:00	06:00
					Epis	sode Start	Time (hh:	mm)		

AT/AF Rx	Rx1:Off	Rx2:Off	Rx3:Off			
Treated episodes						
% Terminated						
ATP Rx versus ms bpm			200-249 241-300			ms bpm
Treated episodes % Terminated						

ATP Sequences: 0 delivered, 0 aborted

### **VT/VF Counters**

Serial Number: PZK875895S

Device: Advisa DR MRI A3DR01 Date of Visit: 07-Mar-2018 13:02:36 Patient: Tsou, Chao Tsung/HA59891 ID: **J120445951** Physician: Liao Min Tsun 0939429152

VT/VF Counters	<b>Prior Session</b> No Prior Session	08-De	Session ec-2017 to ar-2018 lys
VT (>150 bpm)		0	
VT-NS (>4 beats)		0	
Fast A&V		0	
PVC Runs (2-4 beats)		<0.1	per hour
PVC Singles		16.2	per hour
Runs of VRS Paces		0.0	per hour
Single VRS Paces		0.0	per hour

# Arrhythmia Episode List Serial Number: PZK875895S

Device: Advisa DR MRI A3DR01 Date of Visit: 07-Mar-2018 13:02:36 Patient: Tsou, Chao Tsung/HA59891 ID: **J120445951** Physician: Liao Min Tsun 0939429152

Arrhythmia Episode List: 08-Dec-2017 12:46:05 to 07-Mar-2018 13:02:36 All collected episodes.

Type ATP Succe	ess ID# Date	Time hh:mm	Duration hh:mm:ss	Avg bpm A/V	Max bpm A/V	Activity at Onset		
AT/AF	2 14-Jan-20	18 16:25	:06	135/108	207/VP	Active		
AT/AF	1 27-Dec-20	11:57	:01:01	130/77	207/VP	Active		
Last	Programmer Sessi	on 08-Dec-20	)17					
(Data prior to last session has not been interrogated.)								


## Monitored AT/AF Episode #1 Serial Number: PZK875895S

Device: Advisa DR MRI A3DR01 Date of Visit: 07-Mar-2018 13:02:36 ID: **J120445951** Patient: Tsou, Chao Tsung/HA59891 Physician: Liao Min Tsun 0939429152

Туре	ATP Seq	Success	ID#	Date	Time hh:mm	Duration hh:mm:ss	Avg bpm A/V	Max bpm A/V	Activity at Onset		
AT/AF			1	27-Dec-2017	11:57	:01:01	130/77	207/VP	Active		
		• V-V		□ A-A			AT/AF	= 350 ms			
						Detection					Term.
	erval (	ms)				•	-		1.0 min —		<del></del>
120 90	0						•			•	
60		•	•	• • • • • • • • • • • • • • • • • • •	•	• • •	•	• • • • • • • • • • • • • • • • • • •	• • •	• • • • • • • • • • • • • • • • • • • •	•
40	0										
20	00 -			- B	<u> </u>	B B B B	l				8
	+		15	10	· · · · · · · · · · · · · · · · · · ·	· · · · · /	//_		45	10 5	
	-20	-	15	-10	-5	5 (	Time (sec	-20 )	-15	-10 -5	0

## Monitored AT/AF Episode #1 Serial Number: PZK875895S

Device: Advisa DR MRI A3DR01 Date of Visit: 07-Mar-2018 13:02:36 ID: **J120445951** Patient: Tsou, Chao Tsung/HA59891 Episode #1 Chart speed: 25.0 mm/sec

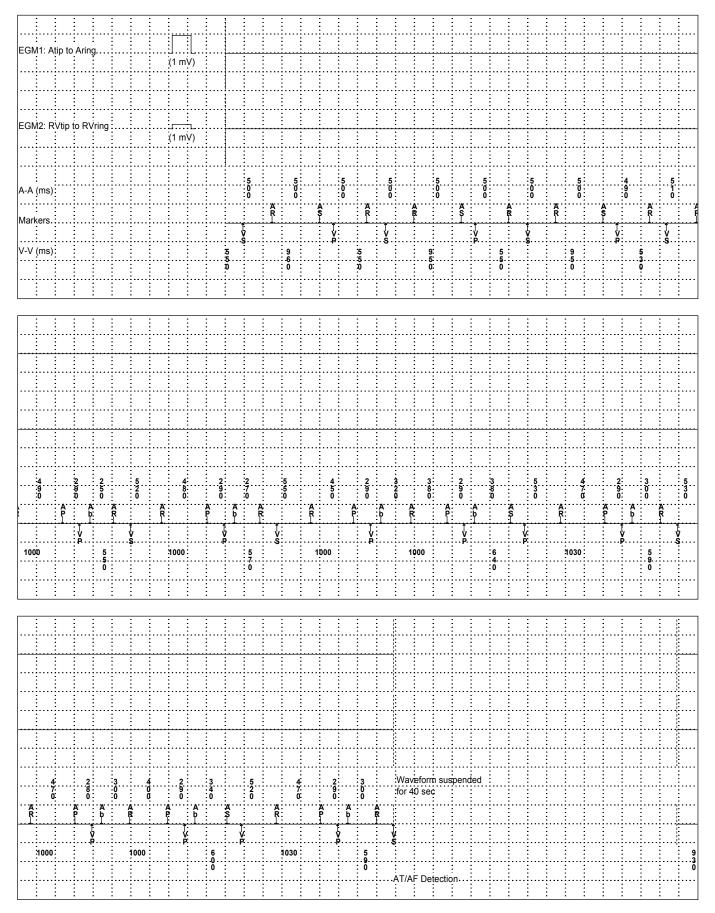


Date of Visit: 07-Mar-2018 13:02:36

Patient: Tsou, Chao Tsung/HA59891

Device: Advisa DR MRI A3DR01

ID: **J120445951** 

Episode #1 Chart speed: 25.0 mm/sec

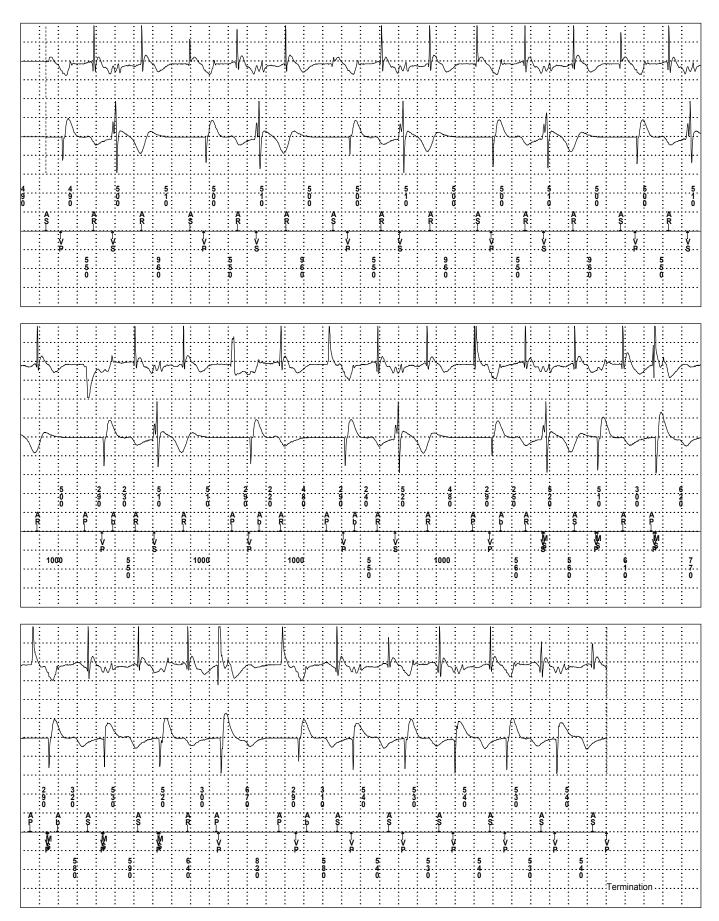

Monitored AT/AF Episode #1 Serial Number: PZK875895S Device: Advisa DR MRI A3DR01 Date of Visit: 07-Mar-2018 13:02:36 Patient: Tsou, Chao Tsung/HA59891 ID: **J120445951** Episode #1 Chart speed: 25.0 mm/sec

						•							•				•					•	-
1 1	: :		: :	: :	: :	:		: :	:		: :	:	:	:	: :			:		:	:	•	:
		:	:	:			:	-:			:	• • • • • • •	• • • • • • • • • • • • • • • • • • • •	:	: :	• • • • • • •		:	: :	• • • • • •	•••••	• • • • • • • • • • • • • • • • • • • •	-:
. i	: :		: :	: :	: :	:	: :	: :	:	•	: :	- 1	:	:	: :			:		:	:		:
	:		:	· · · · · · · · · · · · · · · · · · ·			:	-:				• • • • • • •	• • • • • • • • • • • • • • • • • • • •	:	:	• • • • • • •	• • • • • • • •	:		:	• • • • • • •	• • • • • • • •	:::::
: :	: :	: : :		: :	: :	•		: :	:		: :	:			: :	- 1		:			:	:	•
• • • • • • • • • • • • • • • • • • • •	:	:	:	: : : : : : : : : : : : : : : : : : : :	• : • • • • • • • • • • • • • • • • • •		:	::::::	: : :		: : :	•••••	· · · · · · · · · · · · · · · · · · ·	:	: :	• • • • • •		:		•••••	• • • • • • •	• • • • • • • • •	::
: :	: :	: : :	: :	: :	: :	:	: :	: :	:	:	: :	:	:	:	: :	:	:	: :	: :	:	:	:	:
	: : : : : : : : : : : : : : : : : : : :	: :	:	: :::::::::::::::::::::::::::::::::::::		::::::	: : : : : :	: : : :	:	: : : : : : :	: ::	•••••		:	: ::	: : : :	:	:	:	:	••••		:
: :	: :	: : :	: :	: :	: :	:	: :	: :	:	:	: :	:	:	:	: :	:	:	: :	:	:	:	:	:
: :	: ::::::		:	: : : : : : : : : : : : : : : : : : : :	: : : : :	:::::	: : : : : :	: : :	:	:	: ::	•••••	:	:	: :	•	:	:		••••	••••		:
<del>.</del> <del>.</del>	<b>::</b>		<b>:</b> <del>:</del>	<b>::</b>	. : :	:	: <b>:</b>	.::	<b>:</b>	<b>:</b>	::	: .	<del>:</del>	:	::	<del>:</del>	:	::	: :	:	: .		. <b>:</b>
: :	: :	: :	: :	: :	: :	:	: :	: : :	:	:	: :	:	:	:	: :	:	:	:	:	:	:	:	:
	::		: :	::	. : :	:	: :		: <del>.</del>			: .	<del>.</del>	: :	::			i		:			.:
: :	: :	: : :	: :	: :	: :	:	: :		:		: :	:	:	;	: :			: :	:	:			:
	<u>:</u> :		: :	ii	. : :	:	: : :							· :	::	<del>.</del>		·		:			. <b>:</b>
: :	: :	: : :	: :	: :	: :	:	: :		:		: :		•	:	: :			:	:	:	:		•
	<u>:</u> :		: :	ii	. : :		: <del>.</del>					: .			::			i	:	:			.i
: :	: :		: :	: :	: :	:	: :	: :		:	: :	:	:	:	: :	:	:	:		:	:	:	:
	::		: :	::	.::	:	: :	.::			:		:	:	::			::	: :	:		:	.:
: :	: :	: : :		: :	: :	:		: :	:	•	: :	:	:	:	: :	:	•	:		:	:	•	÷
	·····			·····										·				· · · · ·		• • • • • •		;	٠:٠٠
: :	: :	: : :	: :	: :	: :	:	: :	: :	:		: :	:	:	:	: :	:		:	:	:	:	:	:
	<u>:</u>		<u>:</u> <u>.</u>	<u>:</u>	• • • • • • • • • • • • • • • • • • • •	•••••	: · · · · <u>:</u> · · ·	• • • • • • • •			· • • • • • • • • • • • • • • • • • • •	<del>[</del> .	· · · · <del>[</del> · · · ·		: · · · · <u>:</u>	· · · · <u>:</u> · ·		· · · · ·		• • • • • •		🗄	• • • •
: :	: :	: : :	: :	: :	: :	:	: :	: :	:	:	: :	:	:	:	: :	:	:	:	: :	:	:	:	:
	·····	· · · · · · · · · · · · · · · · · · ·	: · · · · : · · · ·	·····	. : :		: · · · · <del>:</del> · · ·	-::			· <del>!</del> · · · · · <del>!</del>	:		· · · · ·	•	• • • • • • •	:	::		:	···· <del>:</del> ·	:	
: :	÷ ÷	: : :	: :	: :	: :	:	: :	: :	: :	:	; ;		:	:	: :	:	:			:	:	:	:
	·····	· · · · · · · · · · · · · · · · · · ·		·····	• • • • • • • • • • • • • • • • • • • •	:	• • • • • • • • • • • • • • • • • • • •	-::		:	•••••••	• • • • • • •		· : · · · ·	•••••••••••••••••••••••••••••••••••••••	• • • • • • •		:	• • • • •	• • • • • • •			•:••
: :	: :	: : :	: :	: :	: :	:	: :		:		: :		:		: :			: :	: :	:	:	:	:

## Monitored AT/AF Episode #1 Serial Number: PZK875895S

Device: Advisa DR MRI A3DR01 Date of Visit: 07-Mar-2018 13:02:36 ID: **J120445951** Physician: Liao Min Tsun 0939429152 Patient: Tsou, Chao Tsung/HA59891

Episode #1: 27-Dec-2017 11:57:46

### **Episode Summary**

Initial Type AT/AF Monitor (spontaneous)

Duration 1.0 min A/V Max Rate 207 bpm/VP A. Median 176 bpm (340 ms) Activity at onset Active, Sensor = 98 bpm

Parame	eter Settings	Zones		A. Interval (Rate)
AT/AF	Monitor	1	AT/AF	350 ms (171 bpm)

EGM	Source	Range	Sensitiv	/ity
EGM1	Atip to Aring	+/- 2 mV	Atrial	0.3 mV
EGM2	RVtip to RVring	+/- 8 mV	RV	0.9 mV

## Monitored AT/AF Episode #2 Serial Number: PZK875895S

Device: Advisa DR MRI A3DR01 Date of Visit: 07-Mar-2018 13:02:36 ID: **J120445951** Patient: Tsou, Chao Tsung/HA59891 Physician: Liao Min Tsun 0939429152

Туре	ATP Seq	Success	ID#	Date	Time hh:mm	Duration hh:mm:ss	Avg bpm A/V	Max bpm A/V	Activity at Onset		
AT/AF			2	14-Jan-2018	16:25	:06	135/108	207/VP	Active		
		<ul> <li>V-V</li> </ul>		□ A-A			AT/AF	= 350 ms			
								D	etection	Te	erm.
	erval ( 0 ¬	ms)							<del> </del>	6 sec	<b>—</b> ↓
120	0 -										
90		•	•	•	• •	• •	• •	•	• •	_	
60		□•		• _	•	n•	□•			• □	
00		П						•	• • • •	· · · · · · · · · · · · · · · · · ·	<b>D</b>
40	0 -										
	-			00	0				0 0	o o_ o	
20	0 -						C				
	-26	-24	-,	22 -20	-18	-16	-14 - Time (sec	12 -10 c)	-8	-6 -4 -	2 0


Device: Advisa DR MRI A3DR01

Date of Visit: 07-Mar-2018 13:02:36

Patient: Tsou, Chao Tsung/HA59891


ID: **J120445951** 

Episode #2 Chart speed: 25.0 mm/sec

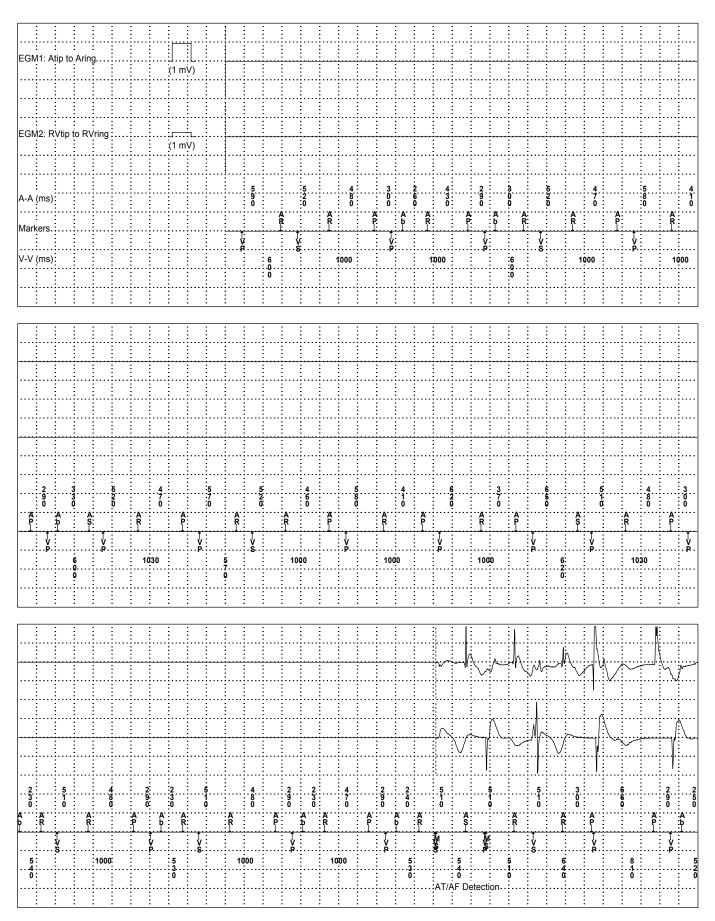

## Monitored AT/AF Episode #2 Serial Number: PZK875895S

Device: Advisa DR MRI A3DR01 Date of Visit: 07-Mar-2018 13:02:36 Patient: Tsou, Chao Tsung/HA59891 ID: **J120445951** Episode #2 Chart speed: 25.0 mm/sec

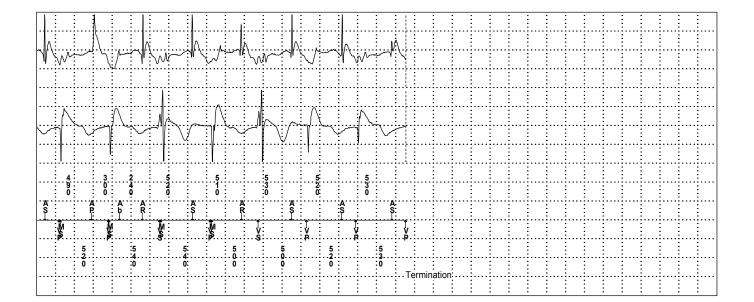

## Monitored AT/AF Episode #2 Serial Number: PZK875895S

Device: Advisa DR MRI A3DR01 Date of Visit: 07-Mar-2018 13:02:36 ID: **J120445951** Physician: Liao Min Tsun 0939429152 Patient: Tsou, Chao Tsung/HA59891

#### Episode #2: 14-Jan-2018 16:25:33

### **Episode Summary**

Initial Type AT/AF Monitor (spontaneous)

Duration 6 sec A/V Max Rate 207 bpm/VP 207 bpm (290 ms) A. Median Activity at onset Active, Sensor = 97 bpm

Parame	eter Settings	Zones		A. Interval (Rate)
AT/AF	Monitor	1	AT/AF	350 ms (171 bpm)

EGM	Source	Range	Sensitivity	
EGM1	Atip to Aring	+/- 2 mV	Atrial	0.3 mV
EGM2	RVtip to RVring	+/- 8 mV	RV	0.9 mV

**Final: Session Summary** 

Device: Advisa DR MRI A3DR01 Serial Number: PZK875895S

Date of Visit: 07-Mar-2018 13:02:36 Physician: Liao Min Tsun 0939429152 Patient: Tsou, Chao Tsung/HA59891 ID: **J120445951** 

**Device Information** 

Device Medtronic Advisa DR MRI A3DR01 PZK875895S Implanted: 08-Dec-2017 Atrial Medtronic 5076 CapsureFix Novus MRI PJN762521G Implanted: 08-Dec-2017 5076 CapsureFix Novus MRI RV Medtronic Implanted: 08-Dec-2017 PJN789201G

Device Status (Implanted: 08-Dec-2017)

Battery Voltage (RRT=2.83V) 3.03 V (07-Mar-2018)

Remaining Longevity 9 years (7.5 - 10 years)

(based on initial interrogation)

Atrial(5076) RV(5076) Lead Impedance 570 ohms 570 ohms

Capture Threshold 0.625 V @ 0.40 ms 0.750 V @ 0.40 ms Measured On 07-Mar-2018 07-Mar-2018 Programmed Amplitude/Pulse Width 1.50 V / 0.40 ms 2.00 V / 0.40 ms

Measured P/R Wave 1.5 mV 12.3 mV **Programmed Sensitivity** 0.30 mV 0.90 mV

**Parameter Summary** 

Mode DDD Lower Rate 60 bpm Paced AV 180 ms Upper Track 130 bpm Sensed AV Mode Switch 171 bpm 150 ms

> Upper Sensor 130 bpm

**Detection** Rates **Therapies** All Rx Off

>171 bpm AT/AF Monitor Monitor >150 bpm

**Changes This Session Session Start Current Value** 

No parameters have been changed during the current session.

Page 1

**Final: Session Summary** Printed: 07-Mar-2018 13:26:53

### **Final: Parameters**

Device: Advisa DR MRI A3DR01 Serial Number: PZK875895S

Date of Visit: 07-Mar-2018 13:02:36 Patient: Tsou, Chao Tsung/HA59891 ID: **J120445951** Physician: Liao Min Tsun 0939429152

**Pacing Summary** 

Pacing Summary					
Mode	Rates		AV Intervals		
Mode DDD	Lower		Paced AV	180 ms	
Mode Switch 171 bpm		I	Sensed AV	150 ms	
	Upper Sensor	130 bpm			
Pacing Details	Atrial	RV			
Amplitude	1.50 V	2.00 V			
Pulse Width	0.40 ms	0.40 ms			
Capture Management	Adaptive	Adaptive			
Amplitude Margin	2.0 X	2.0 X			
Min. Adapted Amplitu	de 1.50 V	2.00 V			
Acute Phase Remain		Off			
Acute Phase Comple	ted 06-Feb-2018	06-Feb-2018			
Sensitivity	0.30 mV	0.90 mV			
Pace Polarity	Bipolar	Bipolar			
Sense Polarity	Bipolar	Bipolar			
Lead Monitor	Adaptive	Adaptive			
Min Limit	200 ohms	200 ohms			
Max Limit	2000 ohms	2000 ohms			
AV Therapies			Arrhythmi	a Interventions	
Rate Adaptive AV	Off		A. Rate Sta		Off
•			A. Preferen		Off
Rate Response			Post Mode	0	Off
ADL Rate	95 bpm		Conducted	AF Response	Off
Optimization	On		V. Rate Sta	bilization	Off
ADL Response	3		A dditional	LEagtures	
Exertion Response	3		Additional		0#
Activity Threshold	Medium Low		Rate Drop I	Response	Off
Activity Acceleration	30 sec		Sleep	Atrial Daoina	Off
Activity Deceleration	Exercise 17		NCAP Inte	Atrial Pacing	On 300 ms
ADL Setpoint UR Setpoint	45		MRI SureS		Off
OK Setpoint	40		PMT Interve		On
Refractory/Blanking			PVC Respo		On
PVARP	Auto		V. Safety P.		On
Minimum PVARP	250 ms		Implant Det	•	Off/Complete
PVAB Interval	150 ms		implant bot		Ciii Compicio
PVAB Method	Partial				
A. Blank Post AP	200 ms				
A. Blank Post AS	100 ms				
V. Blank Post VP	200 ms				
V. Blank Post VS	120 ms				
Device Information					
		a DR MRI A3DR01	PZK875		Implanted: 08-Dec-2017
		CapsureFix Novus MRI			Implanted: 08-Dec-2017
		CapsureFix Novus MRI	PJN789	9201G	Implanted: 08-Dec-2017
Device Configuration II	D: 3-0-0				

Final: Parameters

Printed: 07-Mar-2018 13:27:05

**Final: Parameters** 

Device: Advisa DR MRI A3DR01 Serial Number: PZK875895S

Date of Visit: 07-Mar-2018 13:02:36 Patient: Tsou, Chao Tsung/HA59891 ID: **J120445951** Physician: Liao Min Tsun 0939429152

AT/AF Detection

Detection A. Interval (Rate) Monitor AT/AF No Rx 350 ms 350 ms (171 bpm)

AT/AF Therapies

AT/AF ATP Rx1 Shared A. ATP Rx2 Rx3

AT/AF Rx Status Off Off Off A-A Minimum ATP Interval 150 ms

A. Pacing Amplitude 6 V A. Pacing Pulse Width 1.5 ms

Stop Atrial Therapies After

Disable all A. Rxs if A. lead position is suspect?

**Atrial Lead Position Check** 

No measurement since reset.

**VT Detection** 

Detection V. Interval (Rate)

400 ms (150 bpm) Monitor

**Data Collection Setup** 

Source Range

EGM 1 Atip to Aring +/- 2 mV RVtip to RVring EGM 2 +/- 8 mV EGM 3 Can to RVring +/- 8 mV

Monitored EGM1 and EGM2

Pre-arrhythmia EGM Off

AT/AF Daily Burden 6 hr

Avg. V. Rate During AT/AF 6 hr at 100 bpm

OptiVol Threshold

Device Date/Time 07-Mar-2018 13:04

Holter Telemetry Off